




ACCURA

UNBLOCKER

6F Aspiration-Catheter Set

*This catheter sucks – **with negative pressure to positive results.***

Accura products for a blooming health

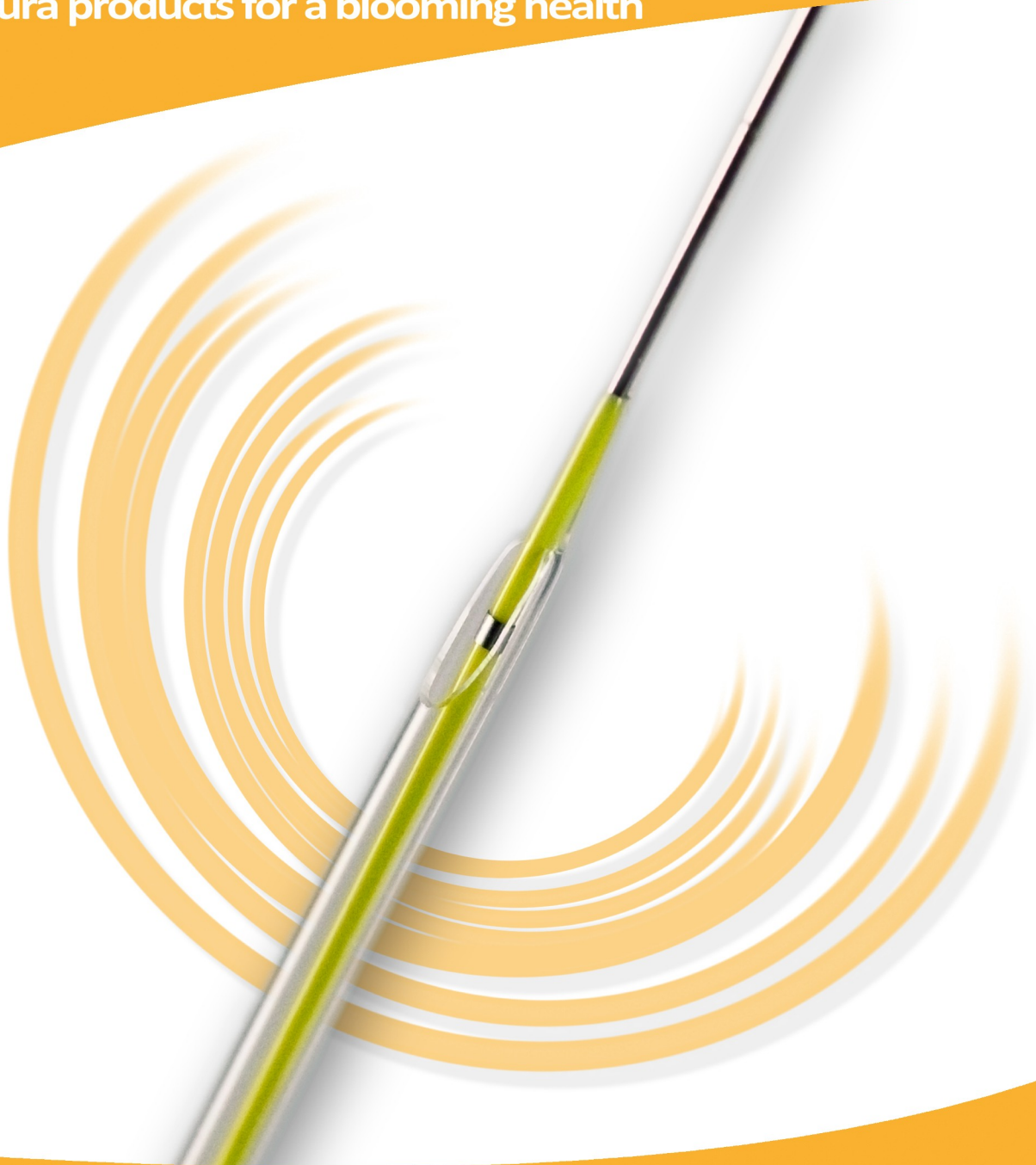

Developed under German Technology
Manufactured in Germany

UNBLOCKER 6F Aspiration Catheter Set


This catheter sucks - *with negative pressure to positive results.*

The UNBLOCKER removes Thrombus and Debris by a low pressure through a **big inner lumen**. The **sophisticated construction** of the catheter gives the catheter a **large inner volumen** while the catheter performs an excellent **pushability** and **tractability**. This all is supported by a **hydrophilic coating** to ease navigation.

Further features are:

- **STD: Short Tip Design** which makes the aspiration procedure easier and more effective
- **Large aspiration port** for an increased aspiration rate
- **Preloaded Stylet** for a better delivery under all conditions and optimal **kink resistance**

It can be applied for **coronary** and **peripheral vessels**.

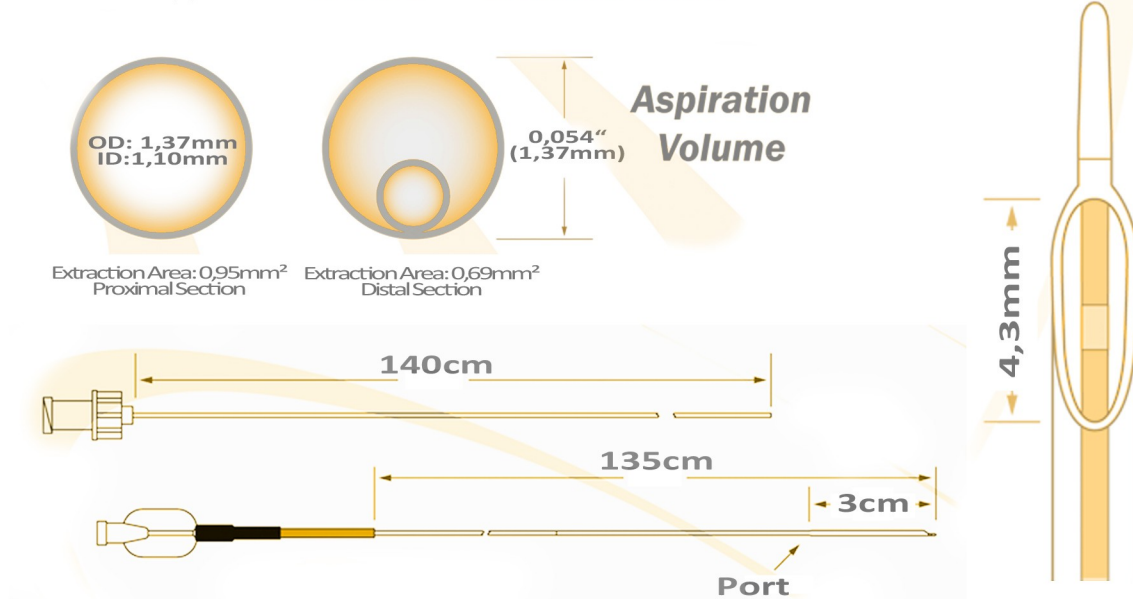

Product Details:

Catheter:

- 135 cm length
- 0.014" Guidewire compatible
- **Rapid Exchange Design**, 30 mm rapid exchange segment
- **Radiopaque marker** 6 mm from tip
- 2 pcs 30 ml **LL Syringes**
- 20 cm **Extension Line with Stopcock**
- 40 μm **Filter Basket**

Catalog Number	Compatible Guide Catheter	Compatible Guide Wire	Shaft Diameter				Guidewire Lumen Part		Effective Length
			Outer		Inner		Max. Outer Diameter	Length	
			Distal	Proximal	Distal	Proximal			
AC6FE	6F	0,014"	4F	4F	1,59F	3,3F	0,572mm	30mm	135cm